


## Halogenated Pesticides Mix # 6 on ZB-50

**Column:** Zebtron™ ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** 50% Phenyl 50% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G004-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min  
**Carrier Gas:** Constant Flow Helium, 0.7 mL/min  
**Injection:** On-Column :1 0.5 µL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes are 10 ppm in Acetone



Products used in this application:

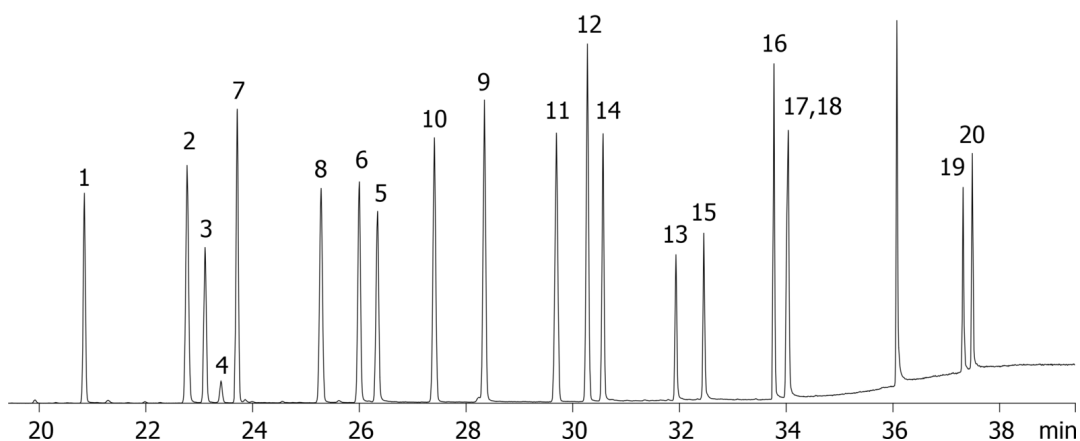

## Halogenated Pesticides Mix # 6 on ZB-50

### ANALYTES:

- 1 Ethalfluralin
- 2 2,3,5,6-Tetrachloroaniline
- 3 Di-allate
- 4 Di-allate
- 5 Simazine
- 6 Atrazine
- 7 Tefluthrin
- 8 Chlordene
- 9 Dimethachlor
- 10 Heptachlor
- 11 Pentachlorothioanisole
- 12 Dacthal
- 13 Pyrifenox
- 14 Fipronil
- 15 Pyrifenox
- 16 Nitrofen
- 17 Etaconazole
- 18 Etaconazole
- 19 Flucythrinate
- 20 Flucythrinate

