

HPLC Application

ID No.: 17884

Terfenadine on Lux 5 μ Cellulose-1 in NP 95% Heptane

Column: Lux[®] 5 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Heptane

Eluent B: 0.1% DEA in IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	10	95	5

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403


- Holder Part No.: KJ0-4282



Products used in this application:



17884

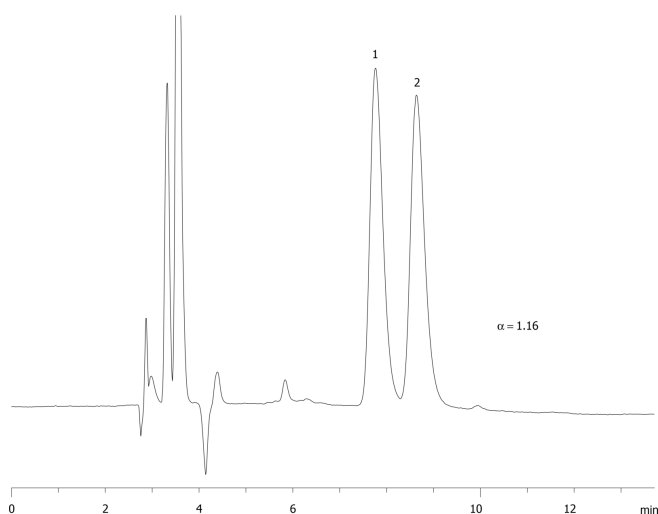

ANALYTES:

- 1 Terfenadine
Retention Time: 7.8 min
- 2 Terfenadine
Retention Time: 8.6 min