

Plasma Metanephries by LC/MS/MS on Luna HILIC following SPE with Strata-X-CW

Column: Luna® 3 µm HILIC 200 Å, LC Column 50 x 2 mm, Ea
Dimensions: 50 x 2 mm ID
Order No: 00B-4449-B0
Elution Type: Gradient
Eluent A: ACN/100mM Ammonium Formate, pH 3.2 (95:5)
Eluent B: ACN/Water/100mM Ammonium Formate, pH 3.2 (50:45:5)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	2.51	0	100
	3	3.01	100	0
	4	6	100	0

Flow Rate: 0.4 mL/min

Col. Temp.: ambient


Detection: Mass Spectrometer (MS) @ amu (ambient)

Detector Info:	Analyte	MRM Pair
		(Q1/Q3)
	MN	180.4/148.2
	MN-D3	183.4/151.2
	NMN	166.4/134.2
	NMN-D3	169.4/137.2

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.
 - SecurityGuard Cartridges, HILIC 4 x 2.0mm, 10/Pk Part No.: AJ0-8328
 - Holder Part No.: KJ0-4282



Products used in this application:

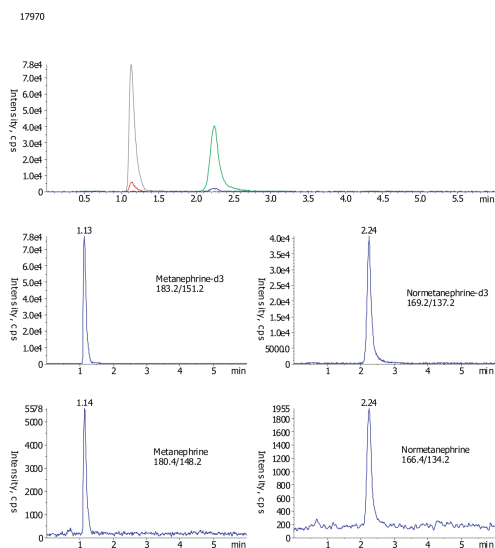

Plasma Metanephrines by LC/MS/MS on Luna HILIC following SPE with Strata-X-CW

ANALYTES:

- 1 Metanephrine
- 2 Metanephrine-d3
- 3 Normetanephrine
- 4 Normetanephrine-d3

