

HPLC Application

ID No.: 18249

Tricyclic Antidepressants on Kinetex 2.6 μ m C18 50 x 4.6 mm ID on Agilent 1100 (system optimization)

Column: Kinetex[®] 2.6 μ m C18 100 \AA , LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4462-E0

Elution Type: Isocratic

Eluent A: 20 mM Potassium Phosphate pH 7.0

Eluent B: 50:50 MeOH and ACN

Gradient

| Step No. | Time (min) | Pct A | Pct B |
|----------|------------|-------|-------|
| 1 | 0 | 48 | 52 |

Profile:

Flow Rate: 2.9 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Detector Info: with 0.12mm ID red capillaries and needle seat, std flow cell, detector rate @ 0.1 s
313 bar



Products used in this application:

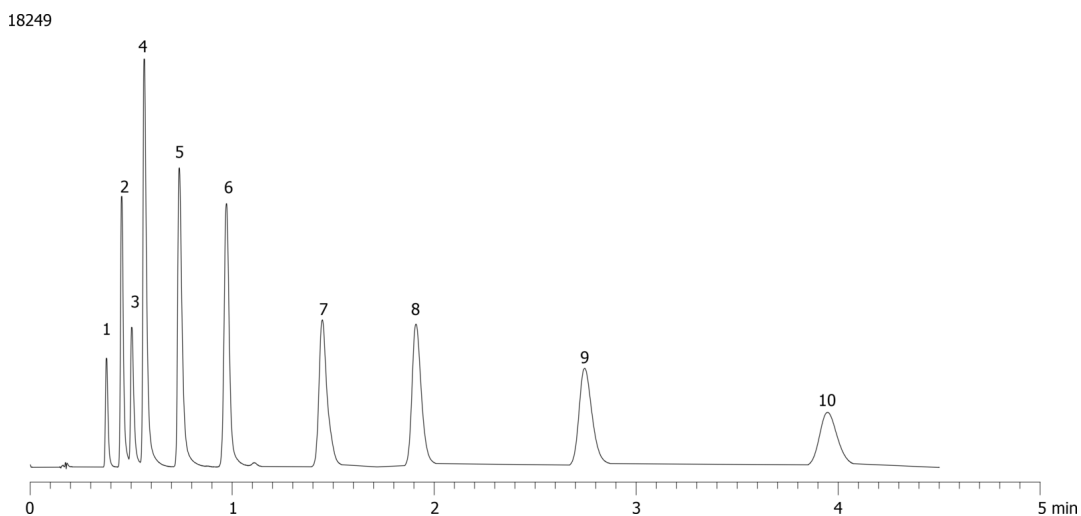

Tricyclic Antidepressants on Kinetex 2.6 μ m C18 50 x 4.6 mm ID on Agilent 1100 (system optimization)

ANALYTES:

- 1 Tianeptine
- 2 Desmethyldoxepin
- 3 Protriptyline
- 4 Desipramine
- 5 Imipramine
- 6 Amoxapine
- 7 Doxepin
- 8 Nortriptyline
- 9 Amitriptyline
- 10 Clomipramine

