

## FMOC-DL-Lysine-(boc)-OH on Lux 5 µm Cellulose-2 250 x 4.6mm in RP

**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4457-E0

**Elution Type:** Isocratic

**Eluent A:** 0.1% TFA in CH<sub>3</sub>CN

**Eluent B:** 0.1% TFA in Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	10	60	40

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282



Products used in this application:

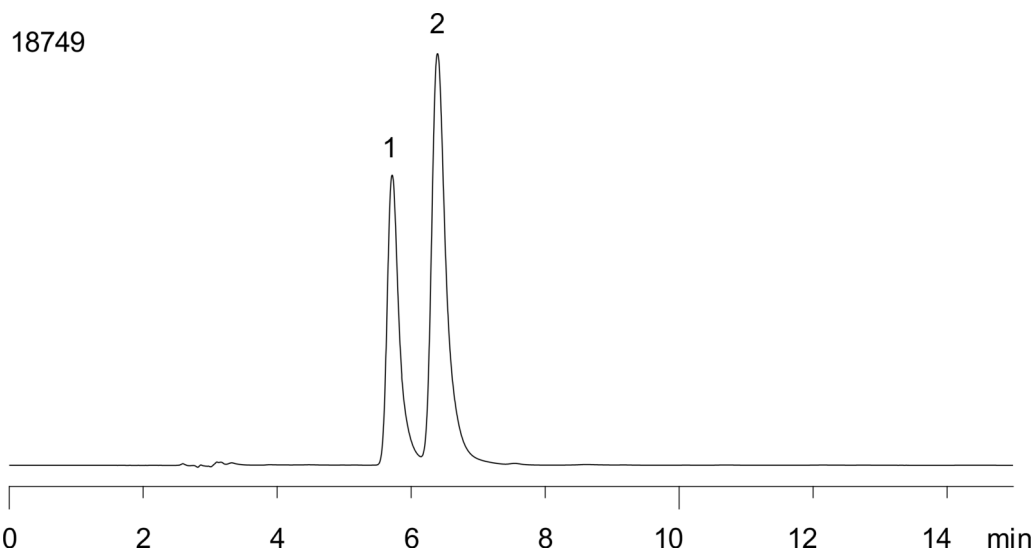

### ANALYTES:

**1** FMOC-Lys(Boc)-OH

Retention Time: 5.7 min

**2** FMOC-Lys(Boc)-OH

Retention Time: 6.3 min

