

Hu IgG2k aggregates on BioSep3000(2) & TSKgel G3000SWxl

Column: BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic

Eluent A: 100mM Sodium Phosphate pH 6.8

Gradient Profile:	Step No.	Time (min)	Pct A
1	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

Detector Info: chromatograms aligned to IgG2 monomer peak


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, GFC 3000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4488

- Holder Part No.: KJ0-4282

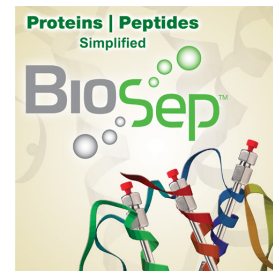

Products used in this application:

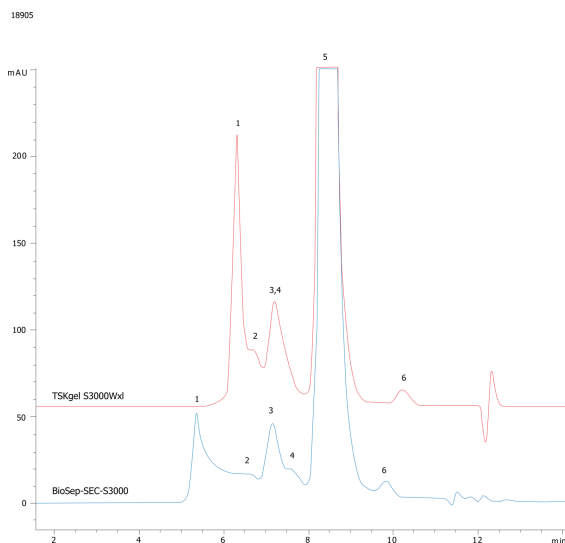

ANALYTES:

- 1** High MW aggregate
Retention Time: 5.2 min
- 2** IgG trimer
- 3** IgG dimer
- 4** IgG2 dimer
Retention Time: 7.1 min
- 5** Human IgG2 kappa monomer
Retention Time: 8.4 min
- 6** impurity

