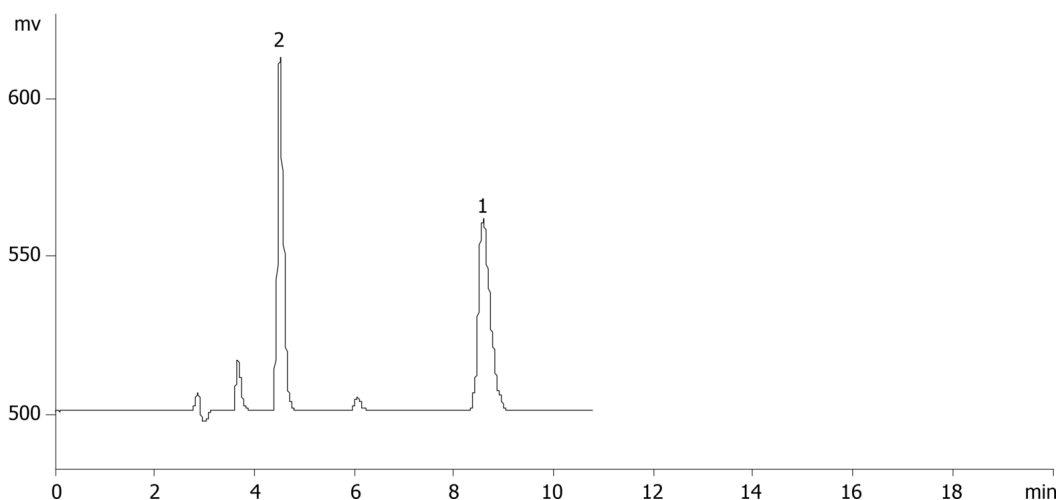


CD Detection of TSO on LUX 5µm Cellulose-2 250 x 4.6 mm ID in NP**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4457-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA in Hexane**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	10	90	10

Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** Circular Dichroism @ 254 nm (25 °C)**Detector Info:** Jasco CDD**Analyst Note:** Please see the corresponding CD trace (attached document).

19150

**Products used in this application:****ANALYTES:**

- 1 trans-Stilbene oxide
- 2 trans-Stilbene oxide

