

## Clarity Oligo-MS: Poly dT 12-18 separation in 5 minutes

**Column:** Clarity® 2.6 µm Oligo-MS 100A, LC Column 50 x 2.1 mm, Ea

**Dimensions:** 50 x 2.1 mm ID

**Order No:** 00B-4479-AN

**Elution Type:** Gradient

**Eluent A:** 4mm TEA/ 200mm HFIP

**Eluent B:** Acetonitrile


| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 100   | 0     |
|                   | 2        | 5          | 70    | 30    |

**Flow Rate:** 0.3 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19174

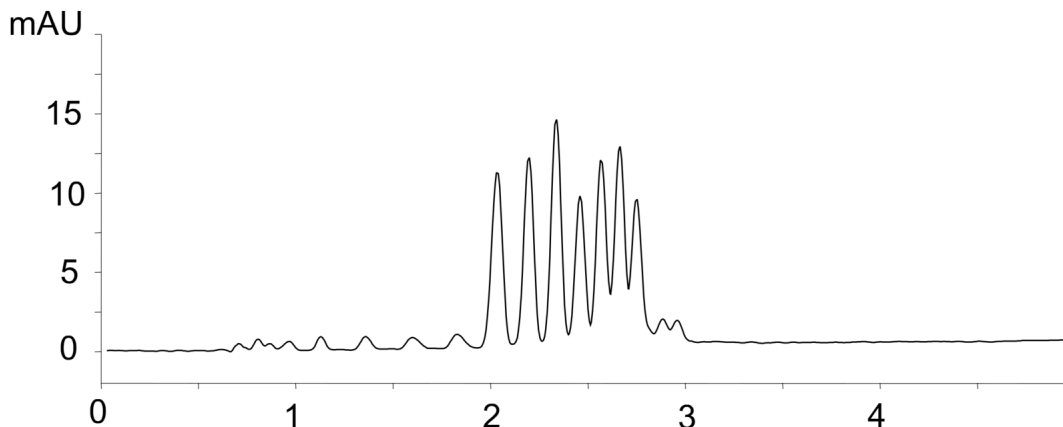

Products used in this application:



### ANALYTES:

- 1 Poly dT 12-18 Oligonucleotide

