

## Poly dT12-18 Oligonucleotides on Luna HILIC

**Column:** Luna® 3 µm HILIC 200 Å, LC Column 50 x 2 mm, Ea

**Dimensions:** 50 x 2 mm ID

**Order No:** 00B-4449-B0

**Elution Type:** Gradient

**Eluent A:** AcCN

**Eluent B:** 10 mM NH<sub>4</sub>OAc

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 90    | 10    |
|                   | 2        | 10         | 50    | 50    |

**Flow Rate:** 0.3 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, HILIC 4 x 2.0mm ID, 10/Pk Part No.: AJ0-8328


- Holder Part No.: KJ0-4282



Products used in this application:



19176

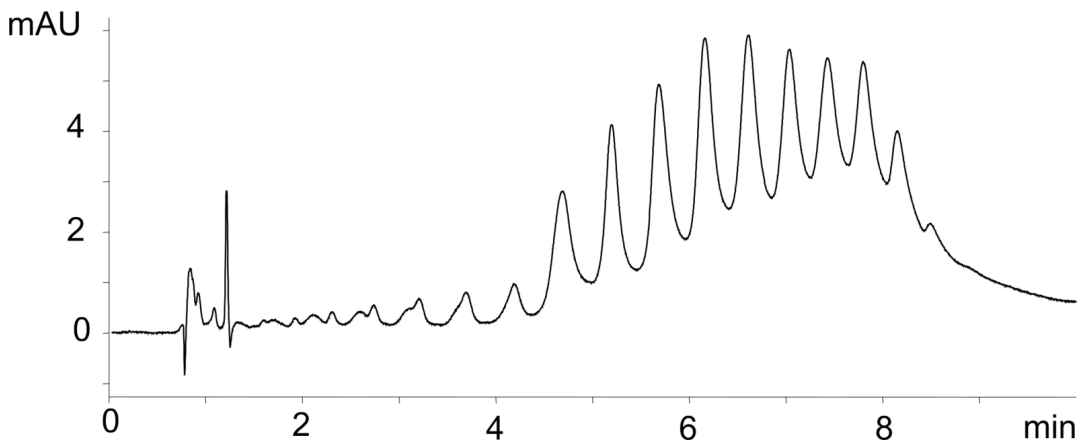

### ANALYTES:

1 poly dT 12-18 mixture

