

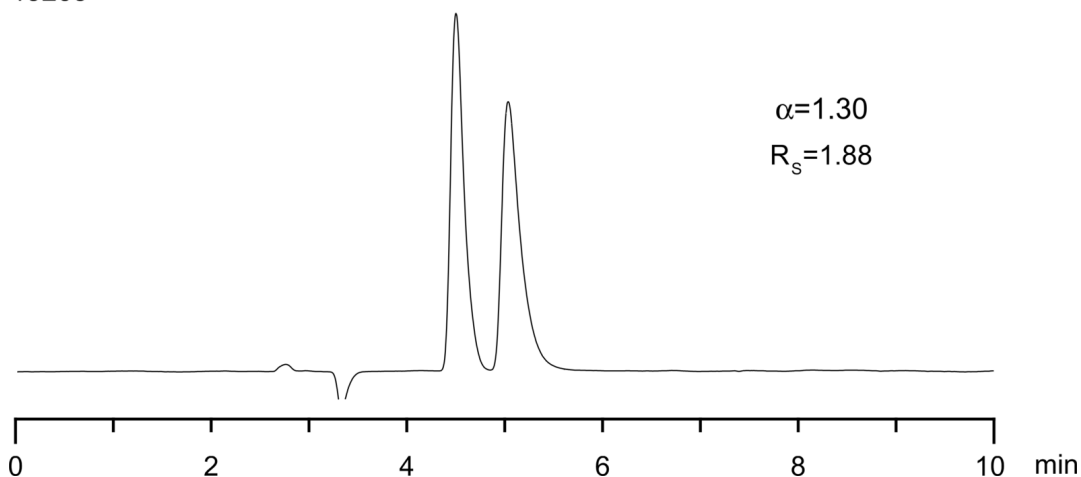
Prilocaine on Lux 5 μ m Cellulose-1 in NP**Column:** Lux[®] 5 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4459-E0**Elution Type:** Isocratic**Eluent A:** Hexane / IPA (80:20) DEA (0.1%)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

**Products used in this application:**

19203

**ANALYTES:****1** Prilocaine

Retention Time: 4.504 min

2 Prilocaine

Retention Time: 5.038 min

