

Chiral HPLC Purification of Terfenadine on Lux Cellulose-2

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 21.2 mm, AXIA™ Packed, Ea

Dimensions: 250 x 21.2 mm ID

Order No: 00G-4457-P0-AX

Elution Type: Isocratic

Eluent A: 0.1% DEA in Isopropanol

Eluent B: 0.1% DEA in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	10	3	97

Flow Rate: 20 mL/min

Col. Temp.: 20 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (20 °C)

Detector Info: PDR Chiral ALP

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard PREP Cartridge, Lux Cellulose-2, 15 x 21.2mm ID, Ea Part No.: AJ0-8400

- Holder Part No.: AJ0-8223

19272



Products used in this application:

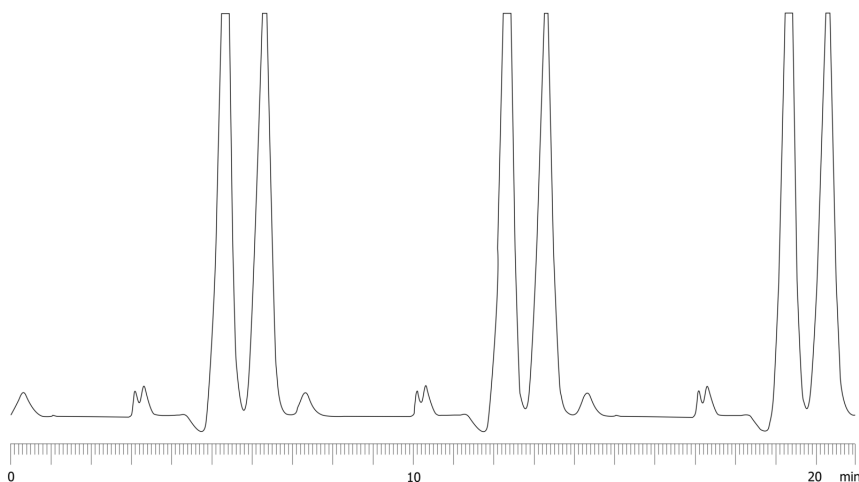

ANALYTES:

1 Terfenadine

Retention Time: 3.5 min

2 Terfenadine

Retention Time: 6 min

