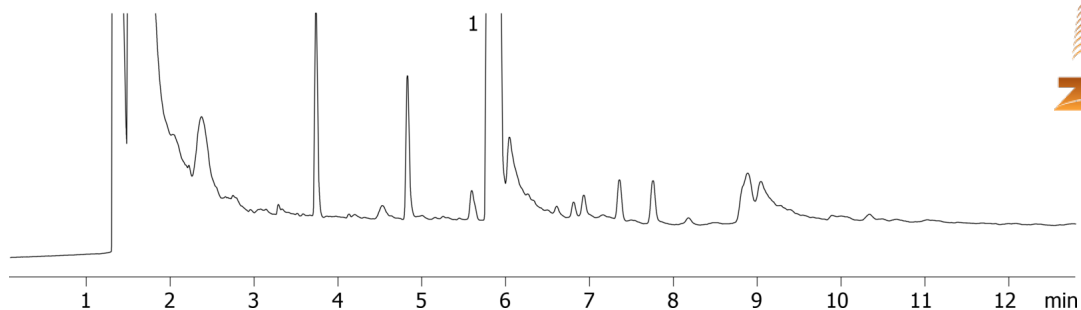


Brompheniramine analysis by FID on a ZB-50 30x0.32x0.5**Column:** Zebtron[™] ZB-50, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G004-17**Oven Profile:** 190°C to 260°C @ 25°C/min for 10 min**Carrier Gas:** Constant Flow Helium, 40 cm/sec**Injection:** Split 20:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (270°C)

19289

**Products used in this application:****ANALYTES:****1** Brompheniramine