


0.1 ng/g Chloramphenicol in Shrimp using Strata-X and Kinetex C18 50 x 2.1 mm by LCMS**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4462-AN**Elution Type:** Gradient**Eluent A:** 5 mM Ammonium bicarbonate**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	2	95	5
	3	2.01	5	95
	4	4.5	5	95

Flow Rate: 0.4 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ (450 °C)**Detector Info:** SCIEX<**Products used in this application:****ANALYTES:**

- 1 Chloramphenicol
Retention Time: 2.1 min
- 2 Chloramphenicol-d5
Retention Time: 2.1 min
- 3 impurity
Retention Time: 3 min

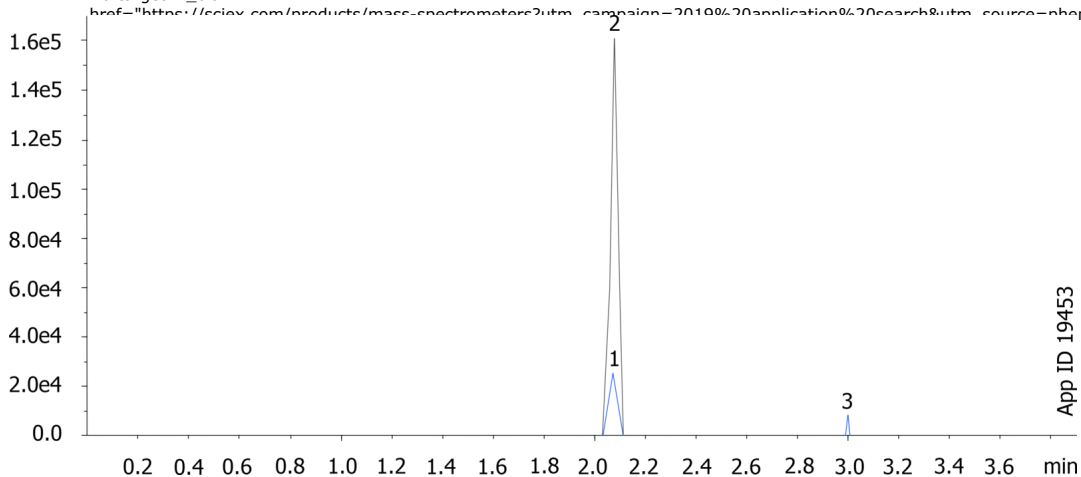