

## USP Loratadine on Luna 5u C8(2) 150x4.60mm

**Column:** Luna® 5 µm C8(2) 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4249-E0

**Elution Type:** Isocratic

**Eluent A:** 31.5:31.5:36.8 Methanol:Acetonitrile:Solution A(1.74mg/mL Anhydrous Dibasic Potassium)

**Gradient**

**Profile:**

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 20         | 100   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282



Products used in this application:

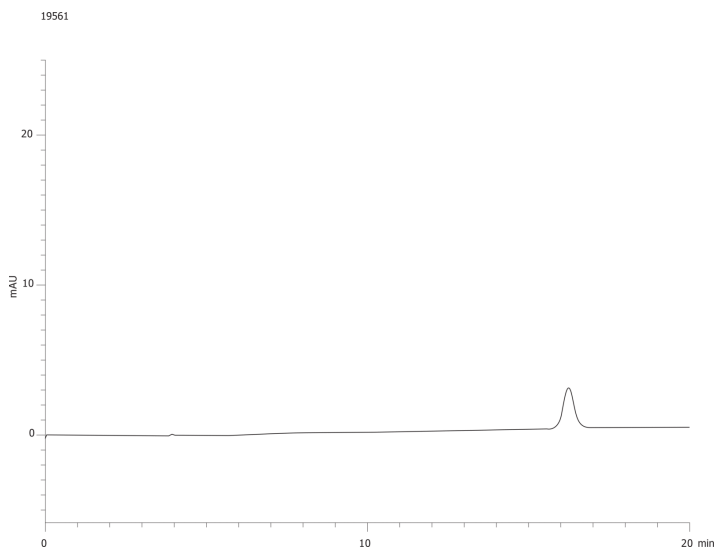

### ANALYTES:

1 Loratadine

