

## Protein Mixture Aeris 3.6µm XB-C8 Wide Pore 150x4.6mm

**Column:** Aeris™ 3.6 µm WIDEPORE XB-C8 200 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4481-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA in Water

**Eluent B:** 0.1% TFA in Acetonitrile

| Gradient        | Step No. | Time (min) | Pct A | Pct B |
|-----------------|----------|------------|-------|-------|
| <b>Profile:</b> | <b>1</b> | 0          | 80    | 20    |
|                 | <b>2</b> | 25         | 35    | 65    |

**Flow Rate:** 1.2 mL/min

**Col. Temp.:** 40 °C

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 214 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C8 for 4.6mm ID, 3/Pk Part No.: AJ0-8771

- Holder Part No.: AJ0-9000



Products used in this application:

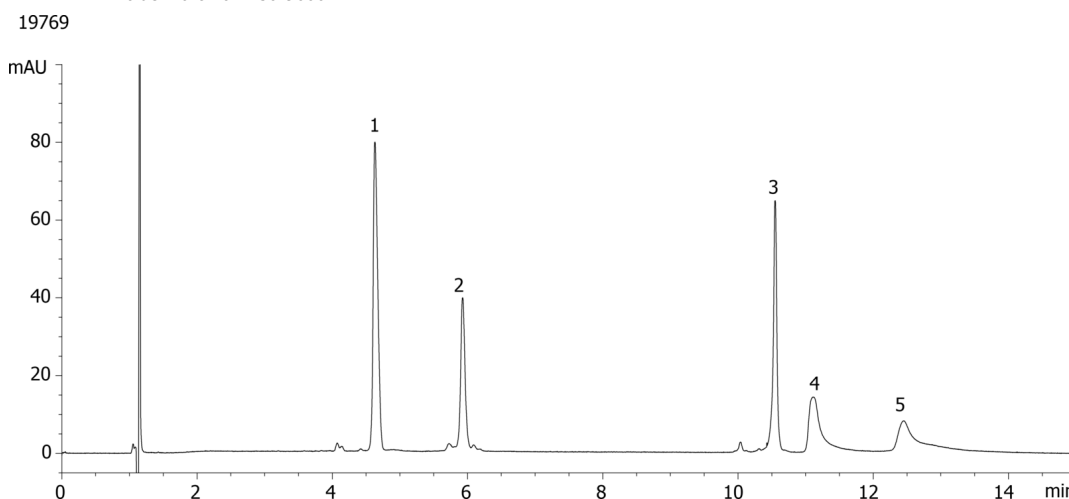

### ANALYTES:

- 1** Angiotensin I  
Retention Time: 4.63 min
- 2** Ribonuclease A  
Retention Time: 5.928 min
- 3** Lysozyme  
Retention Time: 10.55 min
- 4** Apotransferrin  
Retention Time: 11.115 min
- 5** Bovine Serum Albumin (BSA)  
Retention Time: 12.449 min

