

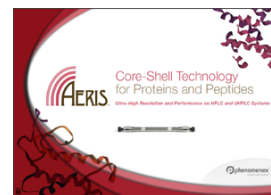
Aeris WIDEPORE 3.6u XB-C18 150 x 2.1mm ID QC Chromatogram (AL0-8931)**Column:** Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 2.1 mm, Ea**Dimensions:** 150 x 2.1 mm ID**Order No:** 00F-4482-AN**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Eluent B:** water

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	55	45

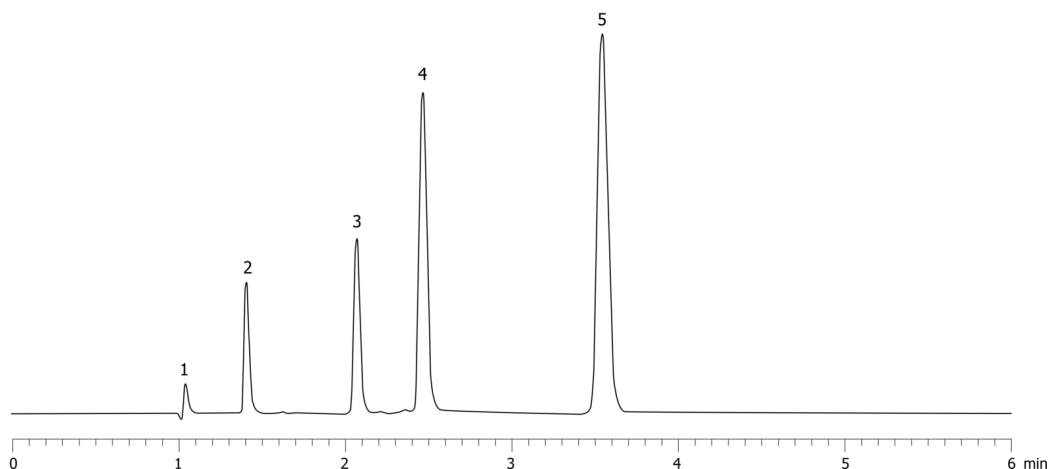
Flow Rate: 0.25 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000

**Products used in this application:**

19812

**ANALYTES:****1** Uracil

Retention Time: 1.041 min

2 Acetophenone

Retention Time: 1.404 min

3 Toluene

Retention Time: 2.066 min

4 Naphthalene

Retention Time: 2.464 min

5 Acenaphthene

Retention Time: 3.543 min

