

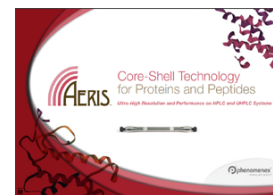
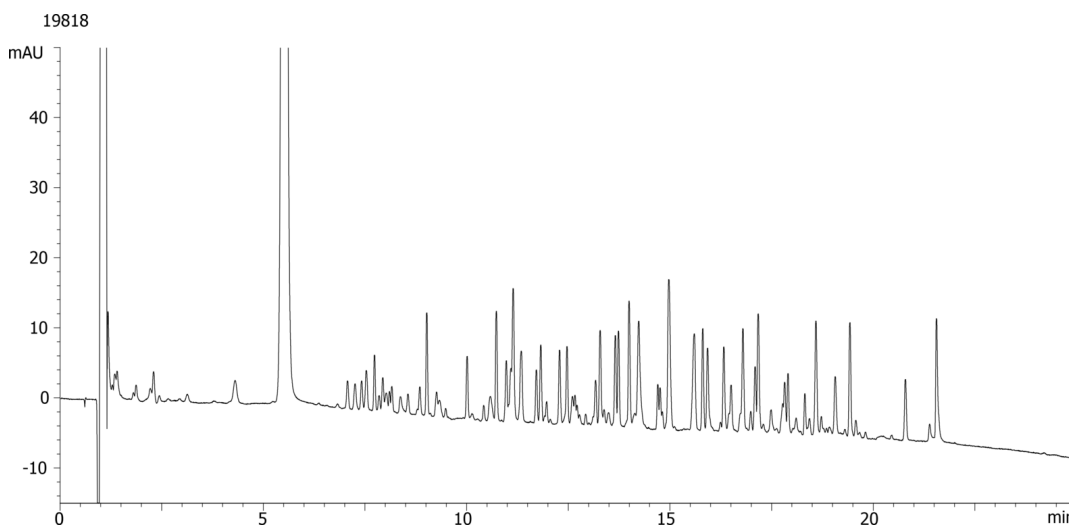
BSA Tryptic Digest on Aeris WIDEPORE 3.6u XB-C18**Column:** Aeris[™] 3.6 μ m WIDEPORE XB-C18 200 Å, LC Column 150 x 2.1 mm, Ea**Dimensions:** 150 x 2.1 mm ID**Order No:** 00F-4482-AN**Elution Type:** Gradient**Eluent A:** 0.1% TFA in Water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35
	2	3	65	35
	3	33	35	65

Flow Rate: 0.3 mL/min**Col. Temp.:** 40 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000

**Products used in this application:****ANALYTES:**

- 1 Bovine Serum Albumin (BSA)

