

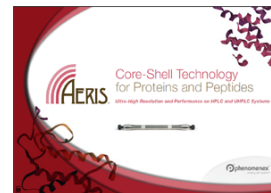
Degraded Myoglobin on Aeris WIDEPORE 3.6u XB-C18 250x2.1 Zoom In**Column:** Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 250 x 2.1 mm, Ea**Dimensions:** 250 x 2.1 mm ID**Order No:** 00G-4482-AN**Elution Type:** Gradient**Eluent A:** 0.1% TFA**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10
	2	5	90	10
	3	50	50	50

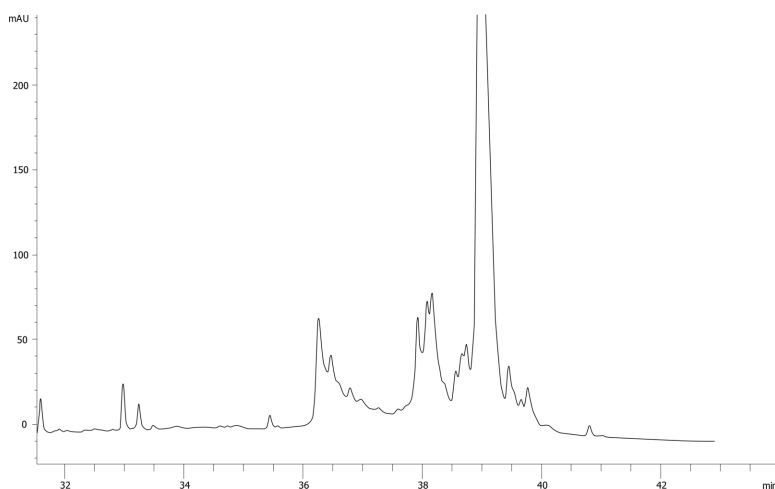
Flow Rate: 0.2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783

- Holder Part No.: AJ0-9000

**Products used in this application:**

19824

**ANALYTES:****1** Myoglobin