

## Apomyoglobin Digest Comparison: Aeris WIDEPORE 3.6µm XB-C18 150x4.6mm (vs19867)

**Column:** Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4482-E0

**Elution Type:** Gradient

**Eluent A:** Water with 0.1% TFA

**Eluent B:** Acetonitrile with 0.1% TFA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	3	97	3
	3	33	35	65

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** ambient


**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 4.6mm ID, 3/Pk Part No.: AJ0-8769

- Holder Part No.: AJ0-9000

19844

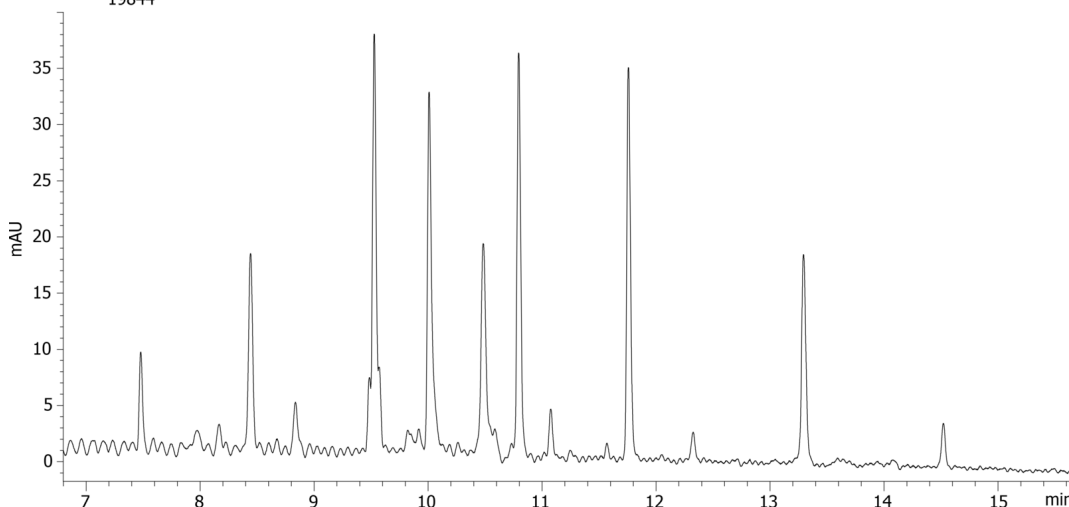

### ANALYTES:

1 Apomyoglobin



Products used in this application:

