

## Security Guard ULTRA - No guard - RNAase A - on Aeris XB-C18 150x4.6mm

**Column:** Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4482-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA/ Water

**Eluent B:** 0.085% TFA/ Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 97    | 3     |
|                   | 2        | 3          | 97    | 3     |
|                   | 3        | 33         | 35    | 65    |

**Flow Rate:** 1.2 mL/min


**Col. Temp.:** 40 °C

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 214 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 4.6mm ID, 3/Pk Part No.: AJ0-8769


- Holder Part No.: AJ0-9000

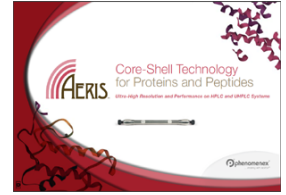

Products used in this application:



19923

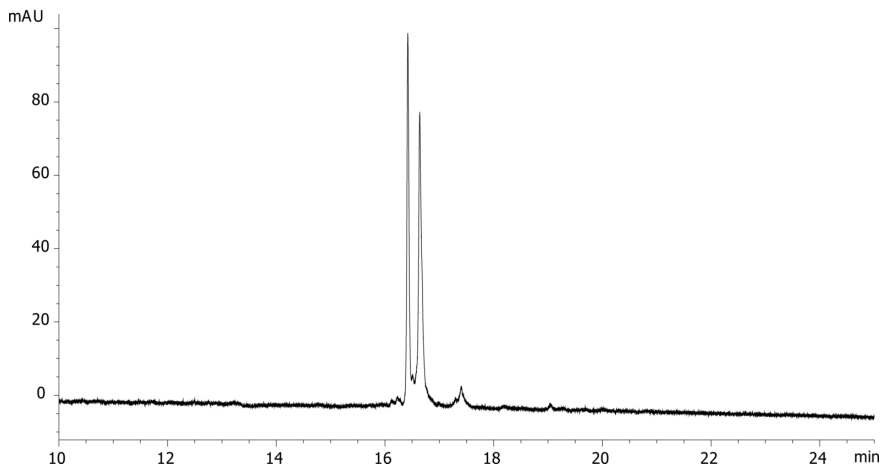

### ANALYTES:

**1** Intact RNase A

Retention Time: 16.42 min

**2** Reduced RNase A

Retention Time: 16.639 min

