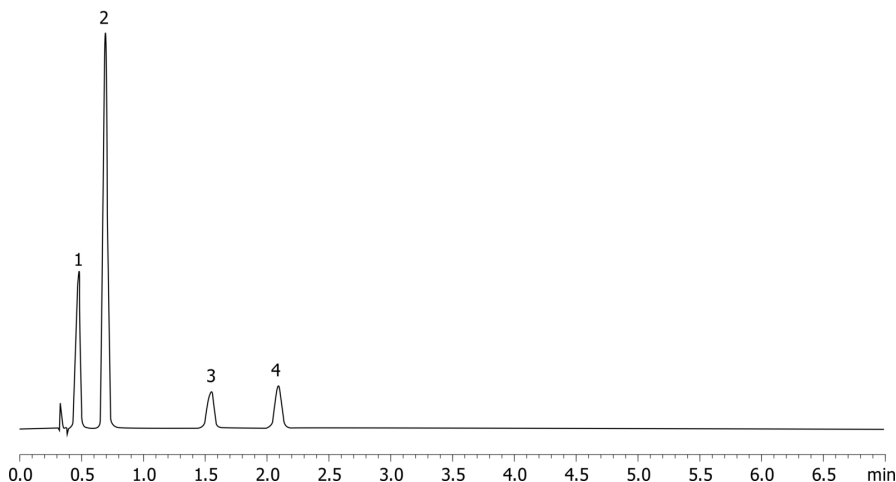


USP Acetaminophen, Acetylsalicylic Acid, Caffeine on Kinetex 2.6u C18 75x4.60mm on Shimadzu Nexera**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 75 x 4.6 mm, Ea**Dimensions:** 75 x 4.6 mm ID**Order No:** 00C-4462-E0**Elution Type:** Isocratic**Eluent A:** 69:28:3 Water:Methanol:Acetic Acid**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 12 100**Flow Rate:** 2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 275 nm (ambient)

19957

**Products used in this application:****ANALYTES:****1** Acetaminophen

Retention Time: 0.488 min

2 Caffeine

Retention Time: 0.706 min

3 Acetylsalicylic acid

Retention Time: 1.564 min

4 Benzoic acid

Retention Time: 2.104 min

