

# HPLC Application

ID No.: 19965

## USP Norethindrone, Progesterone and Mestranol on Kinetex 2.6u C8 150x4.60mm Shimadzu Nexera

**Column:** Kinetex® 2.6 µm C8 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4497-E0

**Elution Type:** Isocratic

**Eluent A:** 50:50 Acetonitrile:Water

**Gradient**

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
|----------|------------|-------|

|                 |          |    |     |
|-----------------|----------|----|-----|
| <b>Profile:</b> | <b>1</b> | 25 | 100 |
|-----------------|----------|----|-----|

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 200 nm (ambient)

19965



Products used in this application:

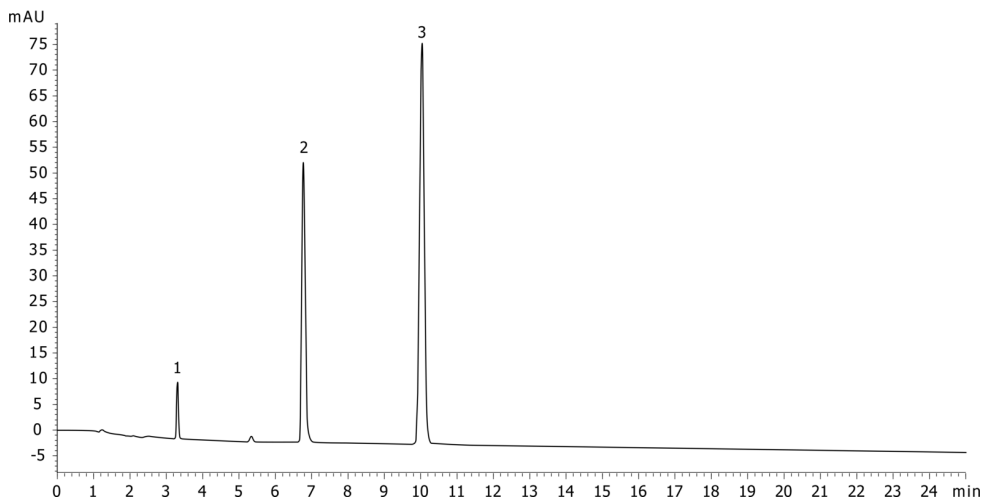

### ANALYTES:

**1** Norethindrone

Retention Time: 3.306 min

**2** Progesterone

Retention Time: 6.775 min

**3** Mestranol

Retention Time: 10.041 min

