

Metomidate on Lux 5µm Cellulose-3 in RP**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4493-E0**Elution Type:** Isocratic**Eluent A:** MeOH with +0.1%DEA**Eluent B:** 20mM NH₄HCO₃

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	80	20

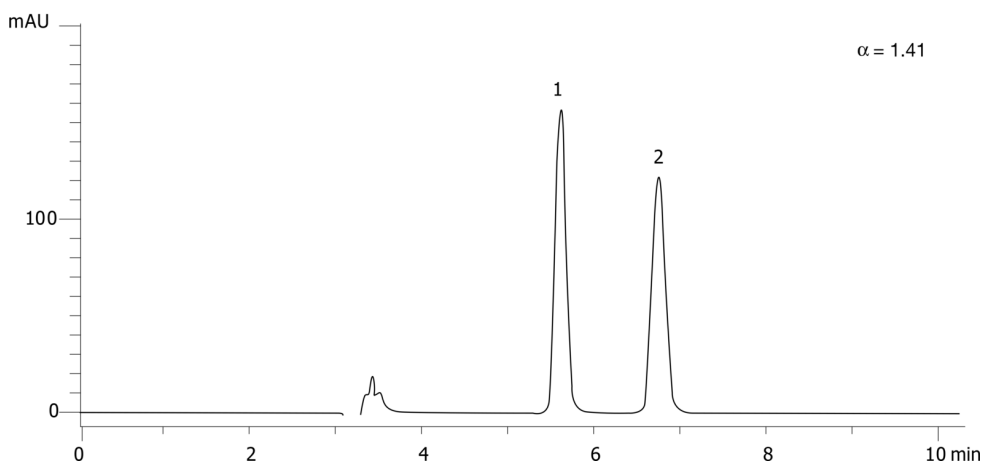
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622

- Holder Part No.: KJ0-4282

**Products used in this application:**

20304

**ANALYTES:****1** Metomidate

Retention Time: 5.62 min

2 Metomidate

Retention Time: 6.75 min

