

Meclizine on Lux 5µm Cellulose-3 in NP**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4493-E0**Elution Type:** Isocratic**Eluent A:** Hexane / EtOH (80:20) DEA (0.1%)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

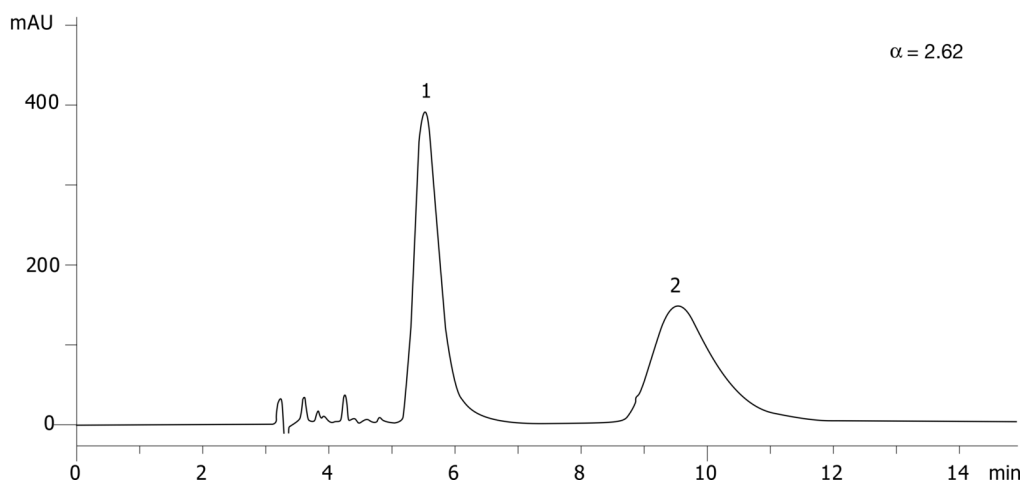
Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622

- Holder Part No.: KJ0-4282

**Products used in this application:**

20338

**ANALYTES:****1** Meclizine

Retention Time: 5.52 min

2 Meclizine

Retention Time: 9.53 min

