

Mirtazapine on Lux 5µm Cellulose-4 in NP**Column:** Lux® 5 µm Cellulose-4, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4491-E0**Elution Type:** Isocratic**Eluent A:** Hexane / IPA (80:20) DEA (0.1%)

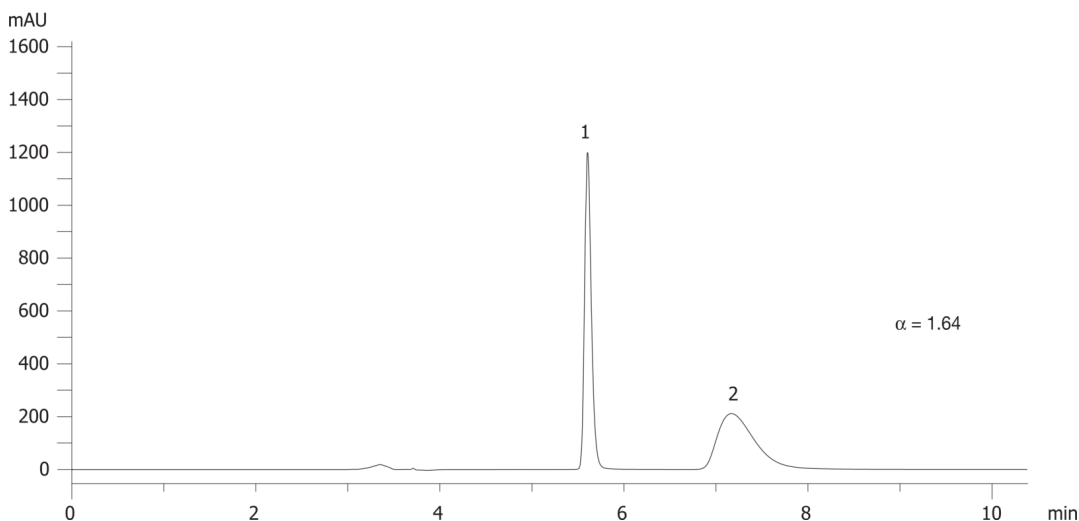
Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-4 4 x 3.0 mm, 10/Pk Part No.: AJ0-8627

- Holder Part No.: KJ0-4282

20425

**Products used in this application:****ANALYTES:****1** Mirtazapine

Retention Time: 5.61 min

2 Mirtazapine

Retention Time: 7.18 min

