

Overlay: Yarra 3µm SEC-3000 and BioSep-SEC-S 3000 Low Salt PMQC Standard

Column: Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-4513-K0

Elution Type: Isocratic


Eluent A: 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% NaN₃

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

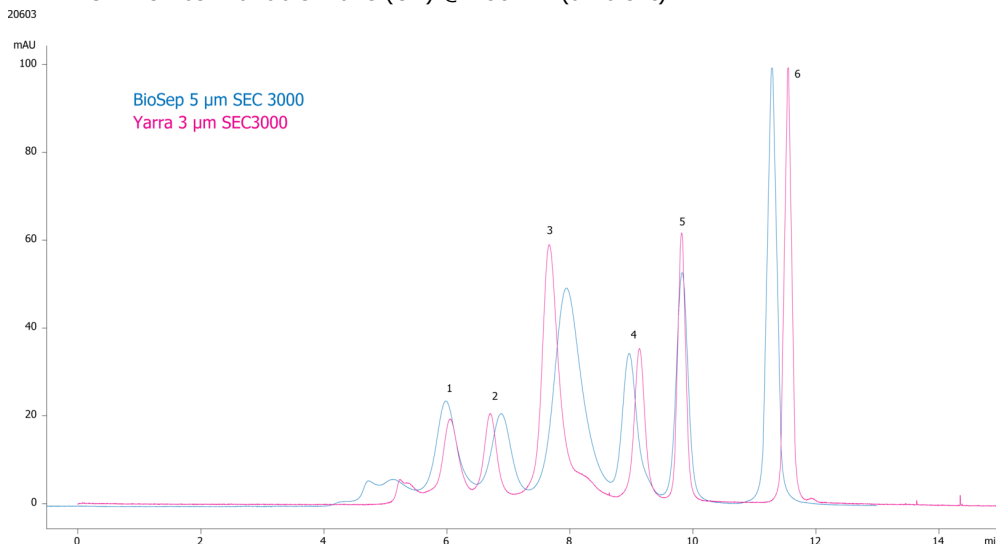

Products used in this application:



ANALYTES:

- | | |
|----------|-------------------------|
| 1 | Thyroglobulin |
| | Retention Time: 6 min |
| 2 | IgA |
| | Retention Time: 6.7 min |
| 3 | IgG |
| | Retention Time: 7.7 min |
| 4 | Ovalbumin |
| | Retention Time: 8.2 min |
| 5 | Myoglobin |
| | Retention Time: 9.8 min |
| 6 | Uridine |
| | Retention Time: 11 min |

