

## Overlay: Yarra 3µm SEC-2000 vs TSKgel super SW2000 low salt buffer PMQC Standard

**Column:** Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-4512-K0

**Elution Type:** Isocratic

**Eluent A:** 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% Na<sub>3</sub>

**Gradient Profile:**


| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

20615

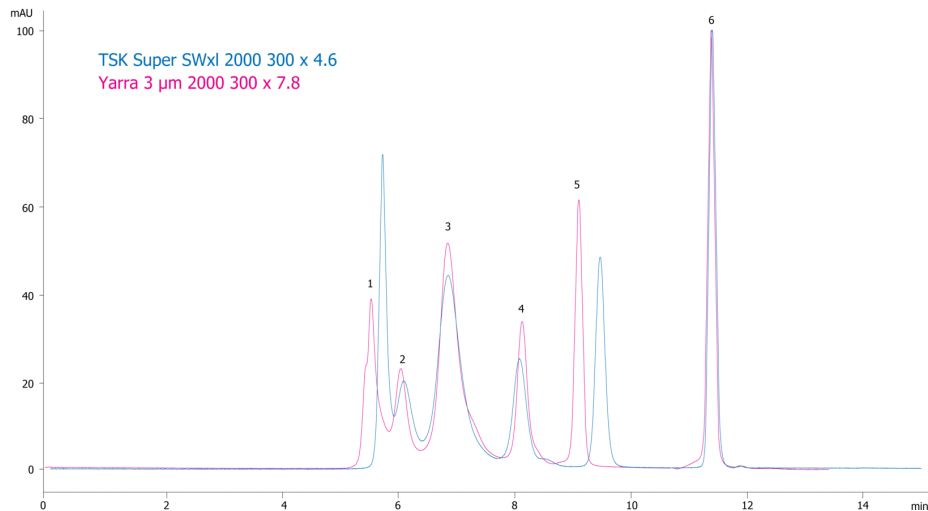

### ANALYTES:

- 1** Thyroglobulin  
Retention Time: 5.5 min
- 2** IgA  
Retention Time: 6.2 min
- 3** IgG  
Retention Time: 7 min
- 4** Ovalbumin  
Retention Time: 8.2 min
- 5** Myoglobin  
Retention Time: 9 min
- 6** Uridine  
Retention Time: 11.5 min



Products used in this application:

