

## Overlay: Yarra 3µm SEC-3000 versus TSKgel 5u G3000 SWxl high salt buffer MMQC Standard


**Column:** Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea  
**Dimensions:** 300 x 7.8 mm ID  
**Order No:** 00H-4513-K0  
**Elution Type:** Isocratic  
**Eluent A:** 50 mM Sodium Phosphate pH 6.8 / 0.3 M Sodium Chloride  
**Gradient**  

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

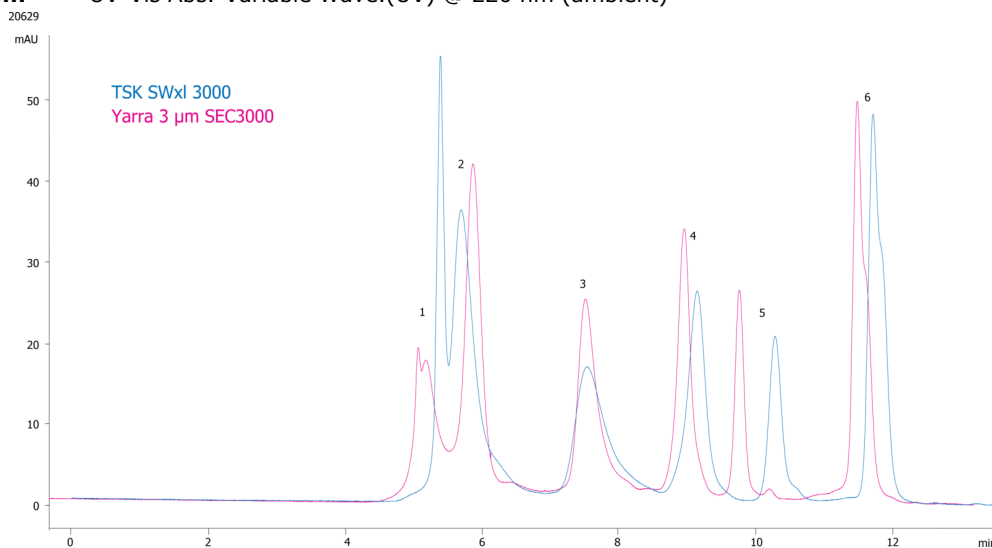

### ANALYTES:

- 1** IgM  
Retention Time: 5.5 min
- 2** Thyroglobulin  
Retention Time: 6.3 min
- 3** Beta Amylase  
Retention Time: 7.9 min
- 4** Ovalbumin  
Retention Time: 9.3 min
- 5** Myoglobin  
Retention Time: 10.1 min
- 6** Uridine  
Retention Time: 11.7 min

