


## Yarra 3µm SEC-2000: G-CSF separations in high salt buffer

**Column:** Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea  
**Dimensions:** 300 x 7.8 mm ID  
**Order No:** 00H-4512-K0  
**Elution Type:** Isocratic  
**Eluent A:** 50 mM Sodium Phosphate pH 6.8 / 0.3 M Sodium Chloride  
**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20641

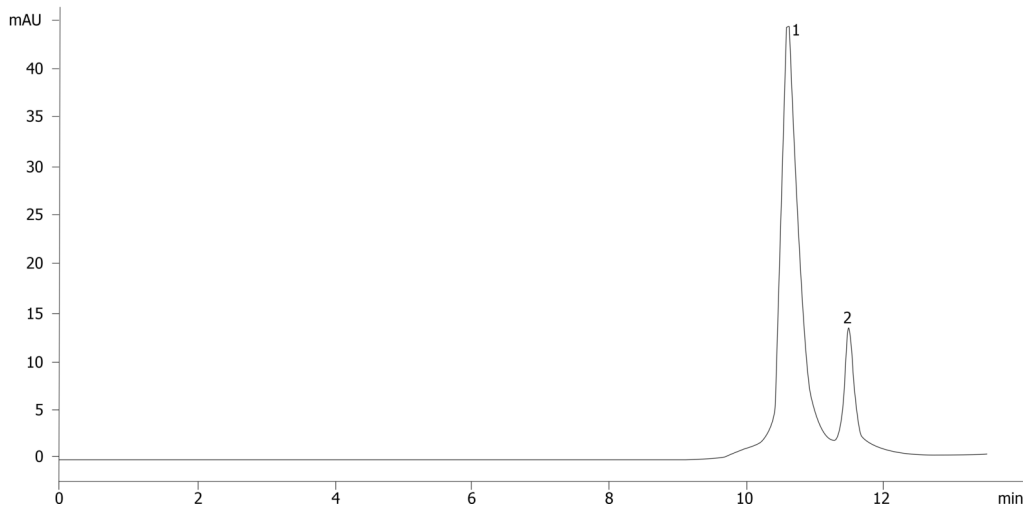

### ANALYTES:

- 1** G-CSF  
Retention Time: 10.49 min
- 2** fragments  
Retention Time: 11.49 min



Products used in this application:

