

## KNX 5 Brochure - Pharma mix analysis using Kinetex 5um C18 150x4.6mm

**Column:** Kinetex® 5µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4601-E0

**Elution Type:** Gradient

**Eluent A:** Water with 0.1% TFA

**Eluent B:** Acetonitrile with 0.1% TFA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 15         | 5     | 95    |

**Flow Rate:** 1.2 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

20738



Products used in this application:

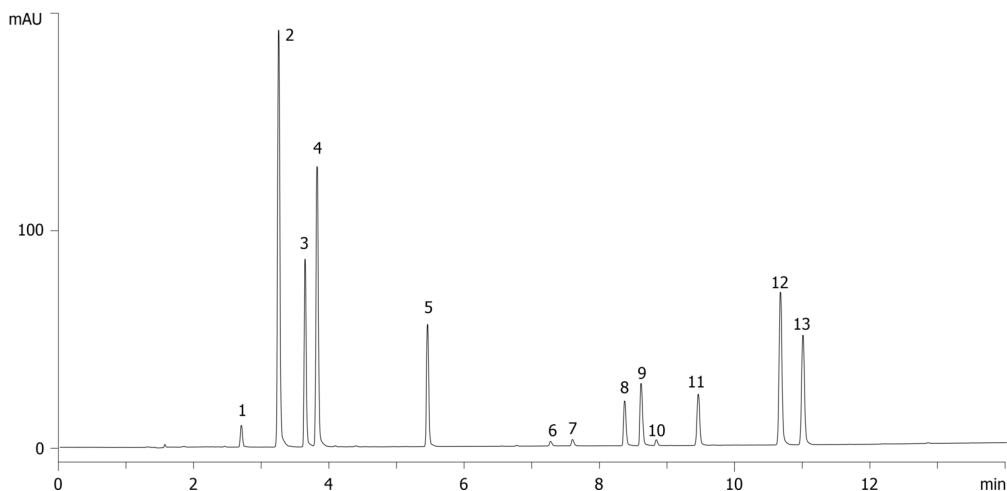

## KNX 5 Brochure - Pharma mix analysis using Kinetex 5um C18 150x4.6mm

### ANALYTES:

- 1** Procainamide
- 2** Acetaminophen
- 3** Folic acid
- 4** Sulfathiazole
- 5** Acebutolol
- 6** Dextromethorphan
- 7** Diphenhydramine
- 8** Propafenone
- 9** Amitriptyline
- 10** Fluoxetine
- 11** Naproxen
- 12** Diflunisal
- 13** Indomethacin

