

LCMS Chiral Separation of Epoxiconazole by Lux 3u Cellulose-2

Column: Lux® 3 µm Cellulose-2, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4456-B0

Elution Type: Isocratic

Eluent A: 5mM Ammonium Hydrogencarbonate

Eluent B: MeOH


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	10	90	

Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: Tandem Mass Spec (MS-MS) @ (ambient)

Detector Info: <a target="_blank"

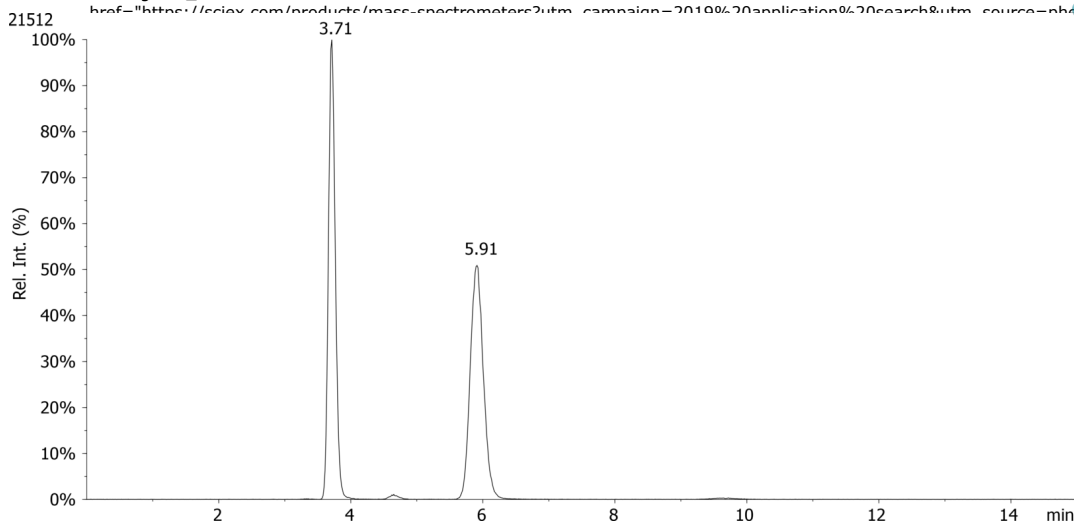

Products used in this application:



ANALYTES:

1 Epoxiconazole

Retention Time: 3.71 min

2 Epoxiconazole

Retention Time: 5.91 min

