

FMOC-Thr(tBu)-OH on Lux 5µm Cellulose-2 in NP

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4457-E0

Elution Type: Isocratic

Eluent A: Hexane / IPA (80:20) Formic acid (0.1%)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

21527



Products used in this application:

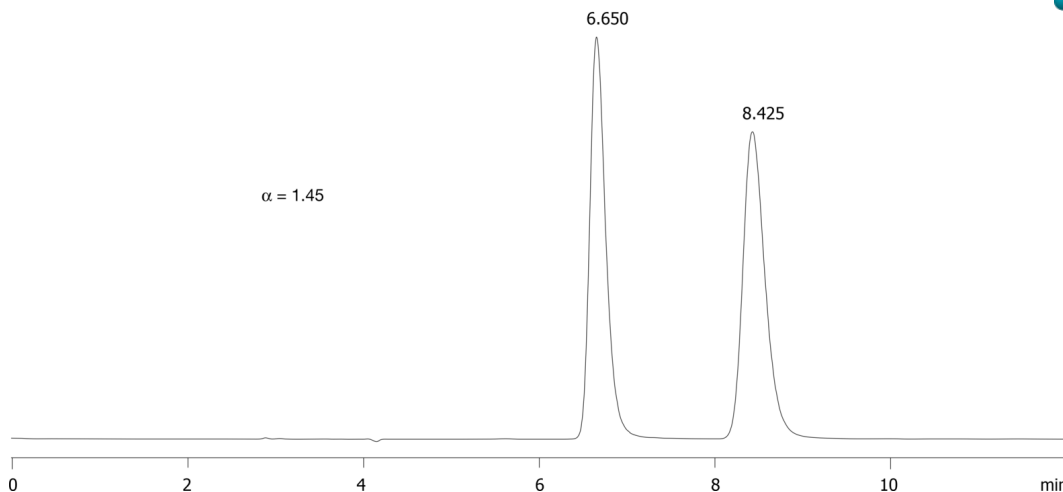

ANALYTES:

- 1** FMOC-Thr(tBu)-OH
Retention Time: 6.65 min
- 2** FMOC-Thr(tBu)-OH
Retention Time: 8.425 min

