

## FMOC-Asp(OtBu)-OH on Lux 5µm Cellulose-3 in RP


**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4493-E0  
**Elution Type:** Isocratic  
**Eluent A:** ACN / 0.1% Formic Acid (40:60) with Formic acid (0.1%) in buffer  
**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)  
21652



Products used in this application:

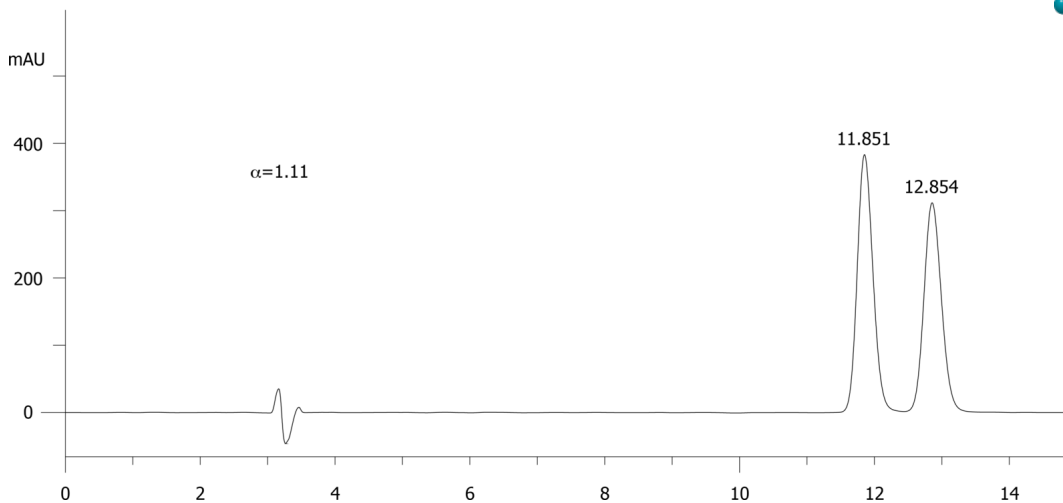

### ANALYTES:

- 1 FMOC-Asp(OtBu)-OH  
Retention Time: 11.851 min
- 2 FMOC-Asp(OtBu)-OH  
Retention Time: 12.854 min

