

# HPLC Application

ID No.: 21792

## Haloxyfop on Lux 5µm Cellulose-3 in RP


**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4493-E0  
**Elution Type:** Isocratic  
**Eluent A:** MeOH / 0.1% Formic Acid (80:20) with Formic acid (0.1%) in buffer  
**Gradient**

Step No.	Time (min)	Pct A
1	0	100

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

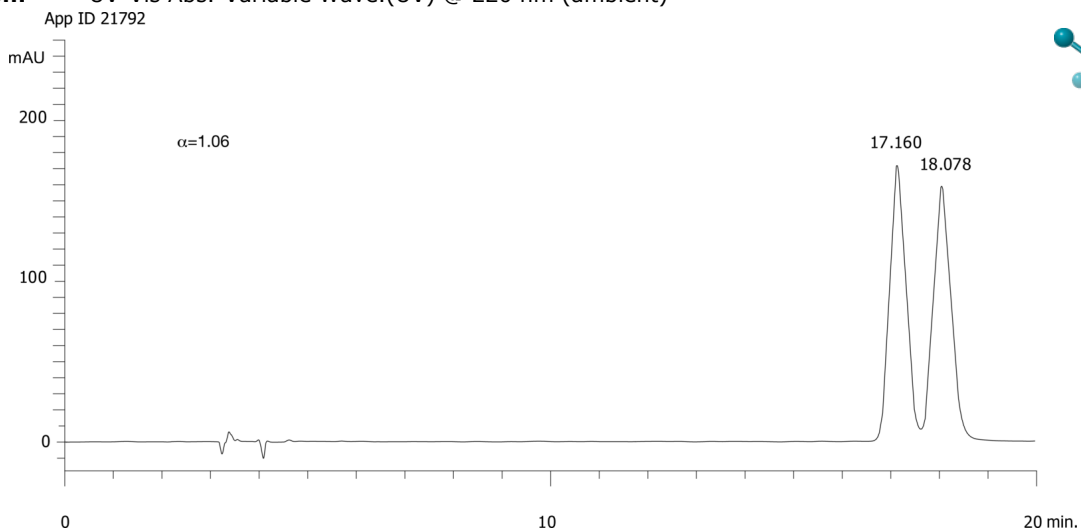

### ANALYTES:

- 1 Haloxyfop  
Retention Time: 17.16 min
- 2 Haloxyfop  
Retention Time: 18.078 min

