# **HPLC Application**

ID No.: 21875



### Human Plasma TMP & Thiamine by Gemini-NX 3u C18 100x3mm with Impact Plate\_2.5nmol

Gemini® 3 µm NX-C18 110 Å, LC Column 100 x 3 mm, Ea

**Dimensions:** 100 x 3 mm ID Order No: 00D-4453-Y0 Elution Type: Gradient

25nM Na2HPO4, 10% methanol Eluent A: Eluent B: 25mM Na2HPO4, 70% methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	3	97
	2	0.25	25	75
	3	0.75	25	75
	4	3	35	65
	5	4	100	0
	6	5	100	0
	7	5.1	3	97
	8	8	3	97




Products used in this application:



Flow Rate: 0.75 mL/min Col. Temp.: ambient

**Detection:** Fluorescence (FLUOR) @ Ex 375 nm (ambient)

**Detector Info:** Shimadzu ELSD

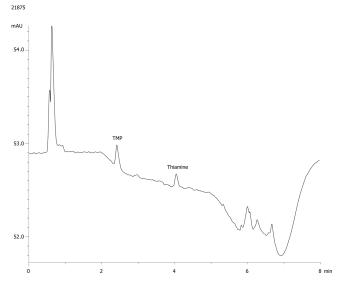

### **ANALYTES:**

TMP

Retention Time: 2.517 min

Thiamine

Retention Time: 4.128 min

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com



## **Sample Preparation Details**

for **HPLC** Application ID No.: 21875



## Human Plasma TMP & Thiamine by Gemini-NX 3u C18 100x3mm with Impact Plate\_2.5nmol

#### **PRODUCT DESCRIPTION:**

Impact™ Protein Precipitation, 2mL Square Well Filter Plate, 2/Pk

Order No.: CE0-7565

### **SOLID PHASE EXTRACTION (SPE) PRODCEDURE:**

**Note:** The solvent volumes shown below are for a Proprietary bed mass.

The solvent volumes will need to be adjusted for a smaller or larger bed mass.

Condition:
Load:
Wash:
Dry:
Elute:
Final Prep and Analysis:
Inject: 20 µL on HPLC Fluorescence (FLUOR) @ Ex 375 nm (ambient)

ANALYTES:	Spiked Conc.	Log P	рКа	% Rec	%RSC
	(ng/mL)				(n=0)
<b>1</b> TMP	2.5				
<b>2</b> Thiamine	2.5				

This method is designed as a convenient starting point for further investigation and can be tailored to meet your extraction goals. Call your local Phenomenex Representative for assistance in method development and optimization techniques.

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com



Phenomenex products are available worldwide. www.phenomenex.com.cn

info@phenomenex.com