

Apo-Myoglobin Tryptic Digest on Aeris PEPTIDE 1.7u XB-C18 150x2.1mm

Column: Aeris™ 1.7 µm PEPTIDE XB-C18 100 Å, LC Column 150 x 2.1 mm, Ea

Dimensions: 150 x 2.1 mm ID

Order No: 00F-4506-AN

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water

Eluent B: 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	1.5	97	3
	3	15	35	65
	4	17	5	95

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22137



Products used in this application:

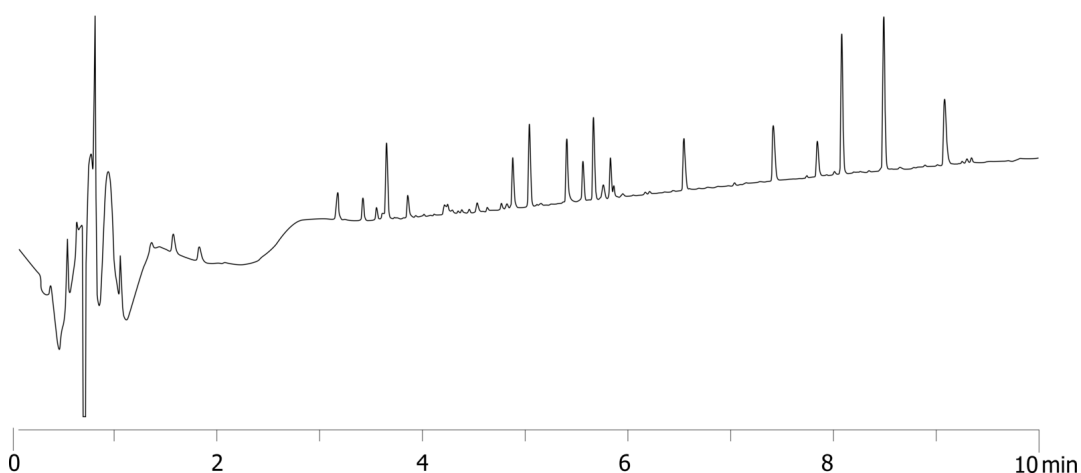

ANALYTES:

- 1 Apo Myoglobin Tryptic Digest

