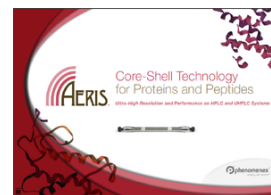
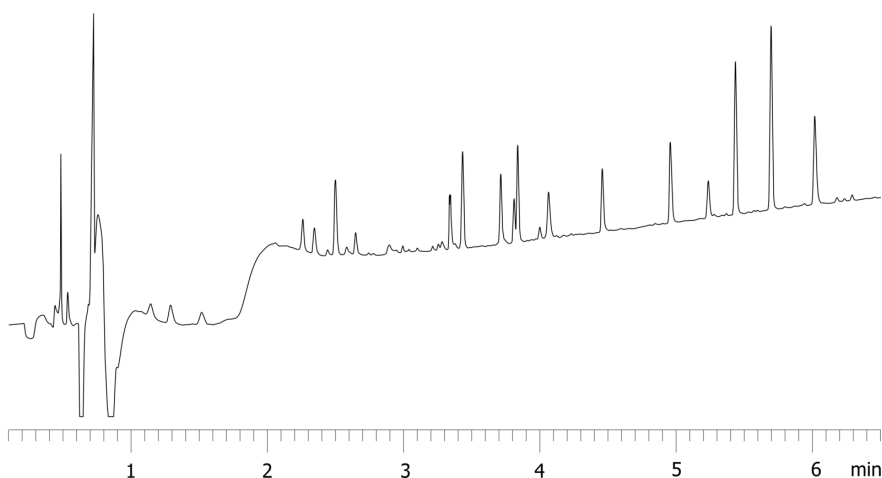


Apo-Myoglobin Tryptic Digest on Aeris PEPTIDE 3.6u XB-C18 250x2.1mm FAST GRADIENT**Column:** Aeris™ 3.6 µm PEPTIDE XB-C18 100 Å, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4507-E0**Elution Type:** Gradient**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	1	97	3
	3	10	35	65
	4	12.5	5	95

Flow Rate: 0.9 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22140

**Products used in this application:****ANALYTES:**

- 1 Apo Myoglobin Tryptic Digest

