

Neonicotinoids using Kinetex 2.6u Biphenyl 50x2.1mm

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4622-AN

Elution Type: Gradient

Eluent A: Water with 0.1% Formic acid

Eluent B: Acetonitrile with 0.1% Formic acid


Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	3	30	70

Flow Rate: 0.6 mL/min

Col. Temp.: 40 °C

Detection: Tandem Mass Spec (MS-MS) @ (ambient)

Detector Info: <a target="_blank"

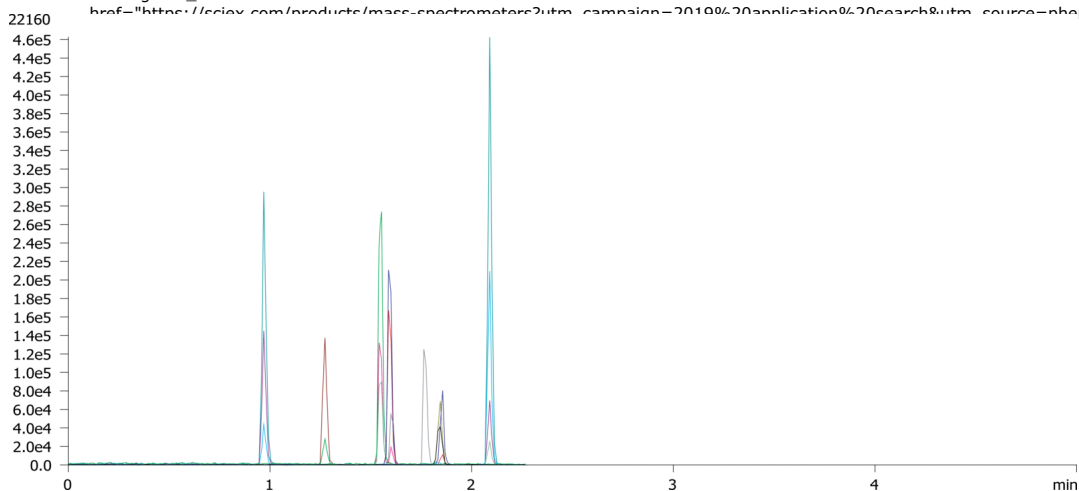

Products used in this application:



ANALYTES:

- 1** Dinotefuran
Retention Time: 0.97 min
- 2** Nitenpyram
Retention Time: 1.27 min
- 3** Thiamethoxam
Retention Time: 1.55 min
- 4** Clothianidin
Retention Time: 1.6 min
- 5** Acetamiprid
Retention Time: 1.85 min
- 6** Thiacloprid
Retention Time: 2.09 min
- 7** Imidacloprid
Retention Time: 2.09 min

