

## TCAs at neutral pH Kinetex 5u EVO C18 150x4.6mm

**Column:** Kinetex® 5µm EVO C18 100 Å, LC Column 150 x 4.6 mm, Ea  
**Dimensions:** 150 x 4.6 mm ID  
**Order No:** 00F-4633-E0  
**Elution Type:** Isocratic  
**Eluent A:** 20mM Sodium phosphate pH 7  
**Eluent B:** Acetonitrile/Methanol 50/50  
**Gradient**


| Step No. | Time (min) | Pct A | Pct B |
|----------|------------|-------|-------|
| 1        | 0          | 30    | 70    |

**Profile:**  
**Flow Rate:** 1.25 mL/min  
**Col. Temp.:** 30 °C  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

22779



Products used in this application:

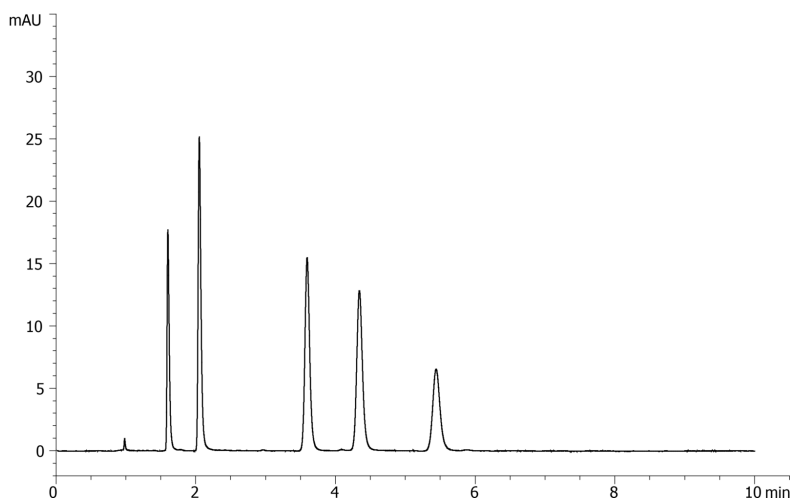

### ANALYTES:

- 1 Protriptyline
- 2 Nortriptyline
- 3 Imipramine
- 4 Amitriptyline
- 5 Clomipramine

