


Pheniramine and Halogenated Derivatives on Kinetex 2.6u EVO C18 100x4.6mm**Column:** Kinetex® 2.6µm EVO C18 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4725-E0**Elution Type:** Gradient**Eluent A:** 10mM Ammonium Bicarbonate, pH=10.2**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50
	2	4	40	60

Flow Rate: 1.25 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 210 nm (25 °C)**Products used in this application:****ANALYTES:**

- 1 Pheniramine
- 2 Chlorpheniramine
- 3 Brompheniramine

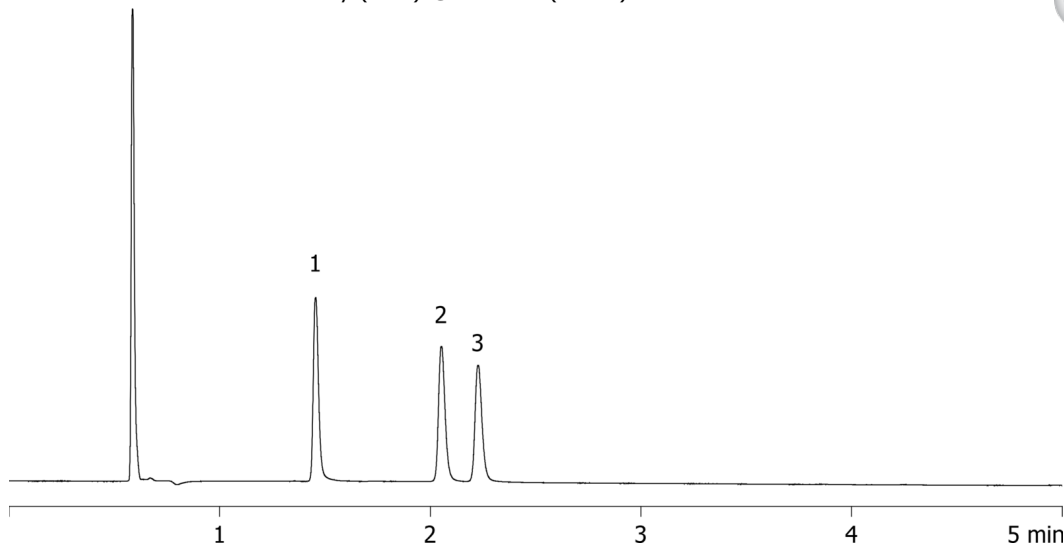