

## Fast Analysis of Sugars in Fruit Juices - Moscato (White Wine)

**Column:** Rezex™ RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

**Dimensions:** 100 x 7.8 mm ID

**Order No:** 00D-0130-K0

**Elution Type:** Isocratic

**Eluent A:** Water

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 0.4 mL/min

**Col. Temp.:** 80 °C


**Detection:** Refractive Index (RI) @ (35 °C)  
23056

**REZEX™**  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX™**

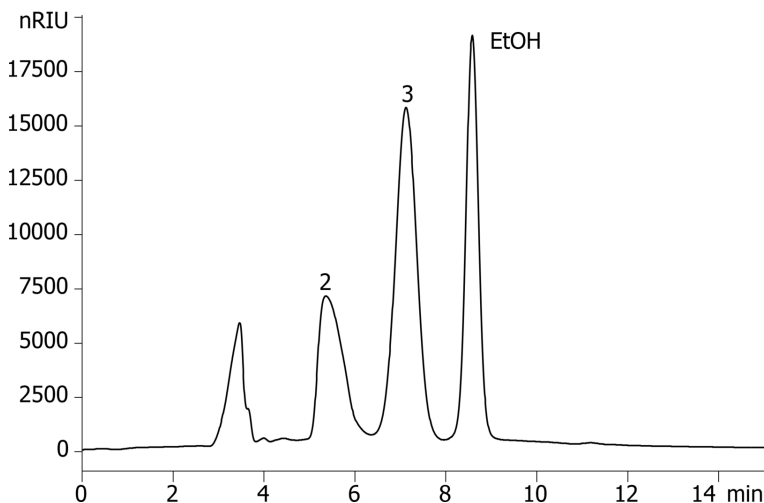

### ANALYTES:

1 Moscato (White Wine)

