


## Yarra 1.8 µm SEC-X150 Batch-to-Batch Reproducibility

**Column:** N/A  
**Dimensions:** N/A  
**Order No:** N/A  
**Elution Type:** N/A  
**Eluent A:** 0.1 M Sodium Phosphate, pH 6.8  
**Gradient**  

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

  
**Profile:**  
**Flow Rate:** 0.3 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

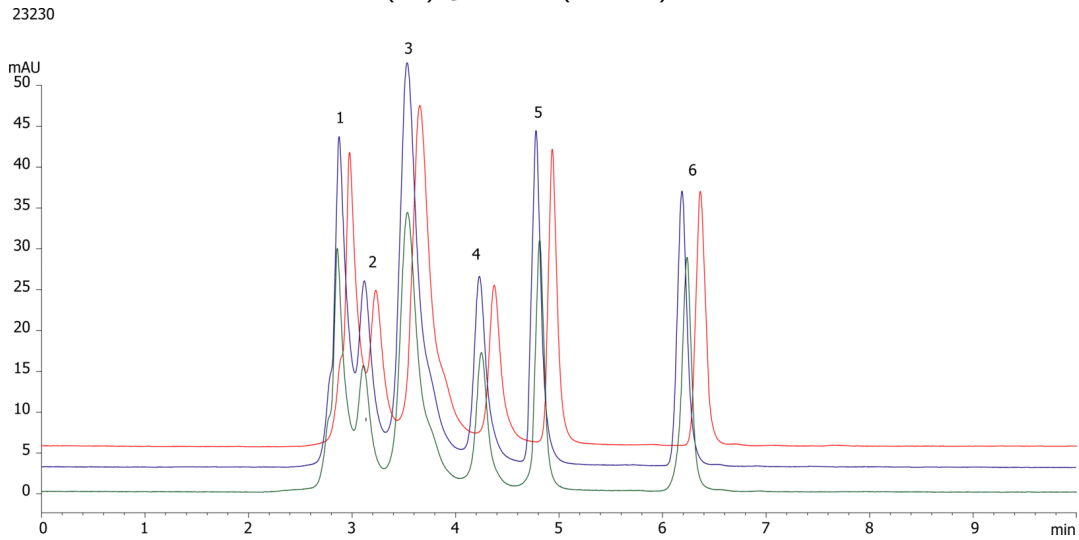

### ANALYTES:

- 1 Thyroglobulin
- 2 IgA
- 3 IgG
- 4 Ovalbumin
- 5 Myoglobin
- 6 Uridine

