


ZB-5MSplus Drugs of Abuse**Column:** Zebron[™] ZB-5MSplus, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G030-11**Oven Profile:** 70 for 1 min to 150 at 30 $^{\circ}$ C/min to 320 at 10 for 5.5 min**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Splitless :1 1 μ L @ 250 $^{\circ}$ C**Detection:** Mass Spectrometer (MS) (250 $^{\circ}$ C)**Detector Info:** 70 for 1 min to 150 at 30 $^{\circ}$ C/min to 320 at 10 for 5.5 min**Analyst Note:** MSD: 45 -450 da
Sample @ 25ppm**Products used in this application:**

23258

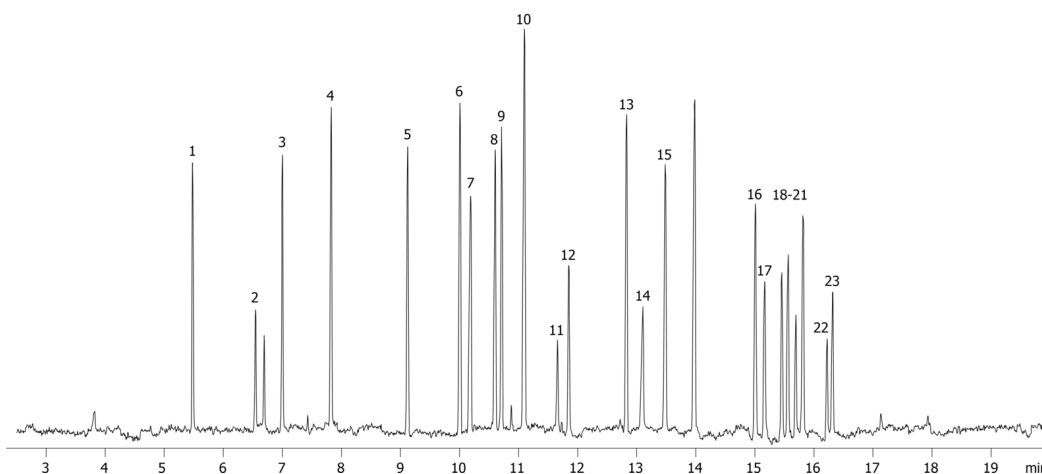

ZB-5MSplus Drugs of Abuse

ANALYTES:

- 1 Nicotine
- 2 MDA
- 3 Mdma
- 4 MBDB
- 5 Meperidine
- 6 Benzamphetamine
- 7 Glutethimide
- 8 PCP
- 9 Doxylamine
- 10 Tramadol
- 11 Procaine
- 12 Methadone
- 13 Topiramate
- 14 Cocaine
- 15 Cocaethylene
- 16 CBD
- 17 Codeine
- 18 Ethylmorphine
- 19 Diazepam
- 20 Hydrocodone
- 21 D9-THC
- 22 6-mam
- 23 Oxycodone

