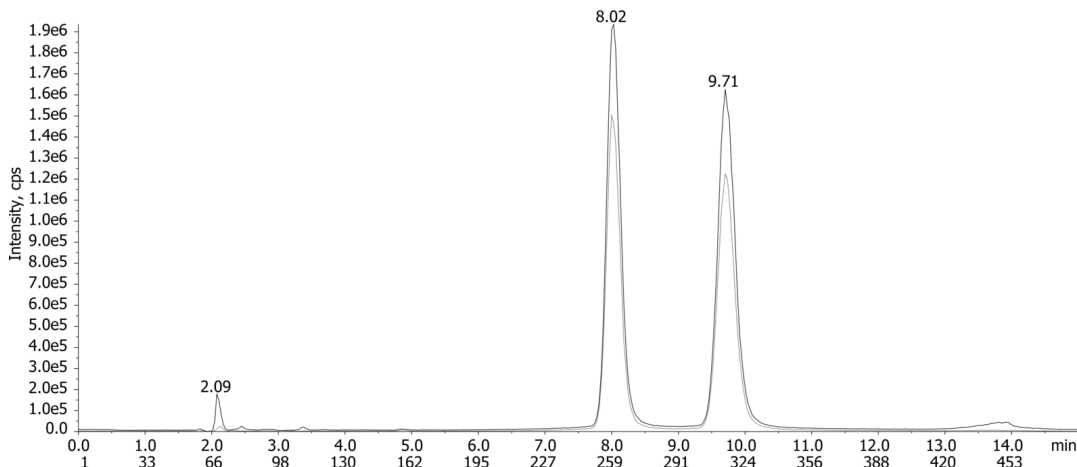


Amphetamine Enantiomers using Lux 3 μ m AMP 150x3mm**Column:** Lux[®] 3 μ m AMP, LC Column 150 x 3.0 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4751-Y0**Elution Type:** Gradient**Eluent A:** 5 mM Ammonium Bicarbonatex**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	40	60
	2	10	40	60
	3	11	5	95
	4	13	5	95
	5	13.1	40	60

Flow Rate: 0.42 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** SCIEX<

23675

**Products used in this application:****ANALYTES:****1** R(-)-Amphetamine

Retention Time: 8.02 min

2 S(+)-Amphetamine

Retention Time: 9.71 min

