

OF Technote Fig 2(ESI-ve): Intercept Device using Knx Biphenyl 50 x 3, 2.6 um

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4622-Y0

Elution Type: Gradient

Eluent A: 10 mM ammonium formate

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50
	2	0.2	50	50
	3	2.5	5	95
	4	4	5	95
	5	4.1	50	50
	6	5	50	50

Flow Rate: 600 µL/min

Col. Temp.: ambient

Detection: Tandem Mass Spec (MS-MS) @ (ambient)

Detector Info: <a target="_blank"


Analyst Note: href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

Analyte panel
Butalbital
Secobarbital
Phenobarbital
THC-COOH

24565



Products used in this application:

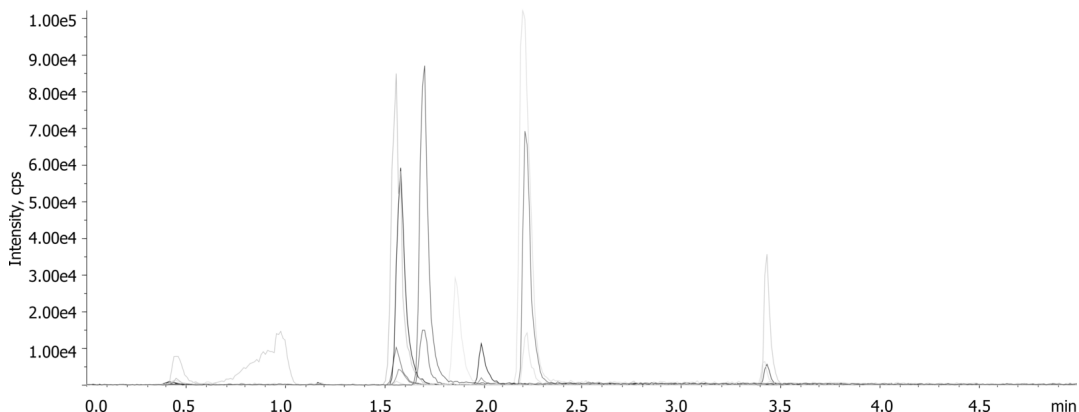

ANALYTES:

1 see analyte panel

