

Analysis of HpCDD, HpCDF, OCDD, OCDF, 37Cl-TCDD and its Isotopes on ZB-Dioxin GC Column
Column: Zebron[™] ZB-Dioxin, GC Cap. Column 60m, 0.25mm, 0.20µm, Ea

Phase:
Dimensions: 60 meters x 0.25 mm x 0.2 µm

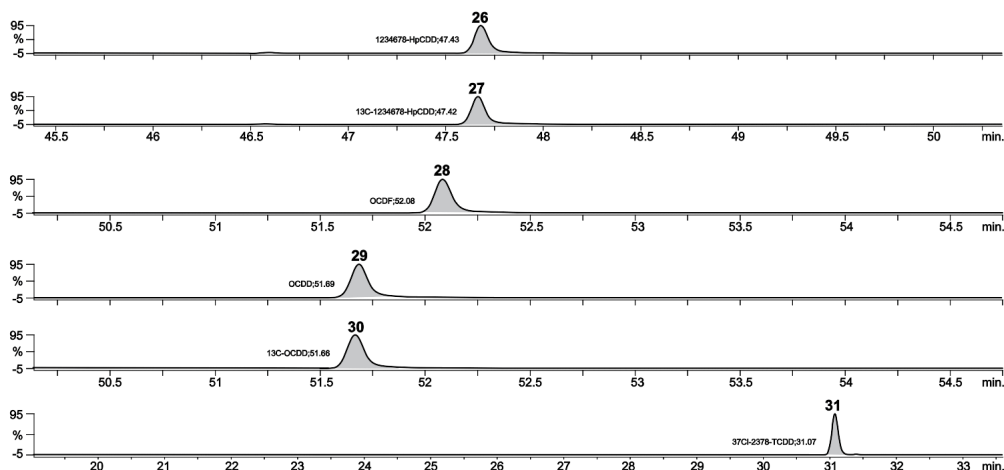
Order No: 7KG-G045-10

Oven Profile: 160°C for 2.4 min to 200°C @ 25°C/min to 220°C @ 5°C/min 220°C for 19 min to 288°C @ 4°C/min to 300°C @ 5°C/min 300°C for 7.6 min

Carrier Gas: Constant Flow Helium, 1.25 mL/min

Injection: Splitless :1 1 µL @ 280°C

Detection: GC/MS (300°C)

Products used in this application:

ANALYTES:

- 1** 1,2,3,4,6,7,8-HpCDD
Retention Time: 47.43 min
- 2** 13C-1,2,3,4,6,7,8-HpCDD
Retention Time: 47.42 min
- 3** OCDF
Retention Time: 52.08 min
- 4** OCDD
Retention Time: 51.69 min

