

PAH with Surrogate Standards on Envirosep-PP (Flu), EPA Method 610

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Water

Eluent B: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	60	40
	2	25	0	100

Flow Rate: 2 mL/min


Col. Temp.: ambient

Detection: Fluorescence (FLUOR) @ 240 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

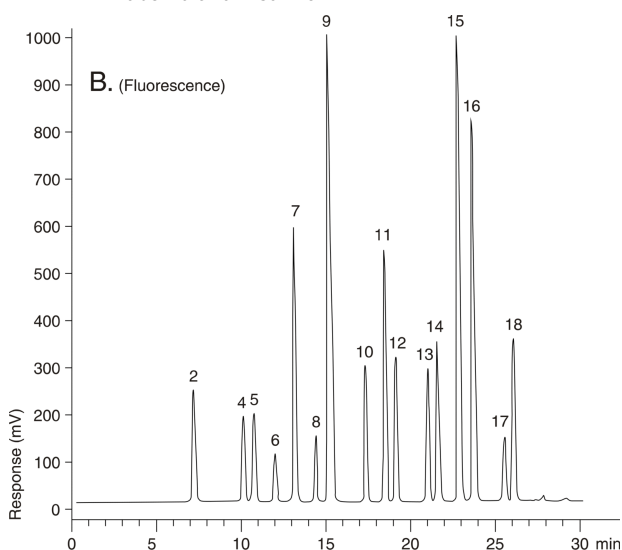

PAH with Surrogate Standards on Envirosep-PP (Flu), EPA Method 610

ANALYTES:

- 1** Nitrobenzene
- 2** Naphthalene
- 3** Acenaphthylene
- 4** Acenaphthene
- 5** Fluorene
- 6** Phenanthrene
- 7** Anthracene
- 8** Fluoranthene
- 9** Pyrene
- 10** p-Terphenyl
- 11** Benzo[a]anthracene
- 12** Chrysene
- 13** Benzo[e]pyrene
- 14** Benzo[b]fluoranthene
- 15** Benzo[k]fluoranthene
- 16** Benzo[a]pyrene
- 17** Dibenzo[a,h]anthracene
- 18** Benzo[g,h,i]perylene
- 19** Indeno[1,2,3-cd]pyrene

