


## Four Cations at 50ppb in the Presence of 1ppm Hydraxine

**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 4 mM Nitric acid in 30% methanol**Gradient**      **Step No.**      **Time (min)**      **Pct A****Profile:**              **1**              0              100**Flow Rate:** 2 mL/min**Col. Temp.:** ambient**Detection:** Conductivity @ nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, RP-1 4 x 3.0mm, 10/Pk Part No.: AJ0-5809

- Holder Part No.: KJ0-4282

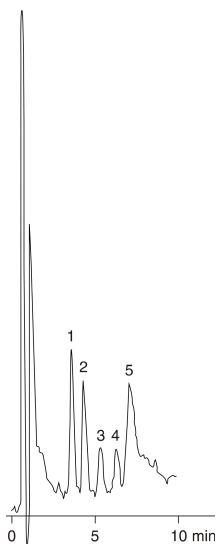

### ANALYTES:

- 1 Lithium
- 2 Sodium
- 3 Ammonium
- 4 Potassium
- 5 Hydrazine

