

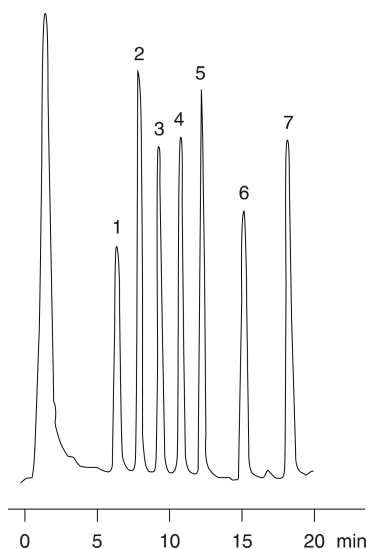
Dansylated Monoamines**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Methanol**Eluent B:** Water mixture

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	55	45
	2	20	90	10

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Fluorimetric @ nm (ambient)**Detector Info:** The following Fluorometric filters were employed: 7-60/7-54 and 3-72/4-76**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Ethylamine
- 2 Propylamine
- 3 Butylamine
- 4 Pentylamine
- 5 Hexylamine
- 6 Octylamine
- 7 Decylamine

