

## Dansylated Monoamines

**Column:** N/A

**Dimensions:** N/A

**Order No:** N/A

**Elution Type:** N/A

**Eluent A:** Methanol

**Eluent B:** Water mixture

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	55	45
	2	20	90	10

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient


**Detection:** Fluorimetric @ nm (ambient)

**Detector Info:** The following Fluorometric filters were employed: 7-60/7-54 and 3-72/4-76

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

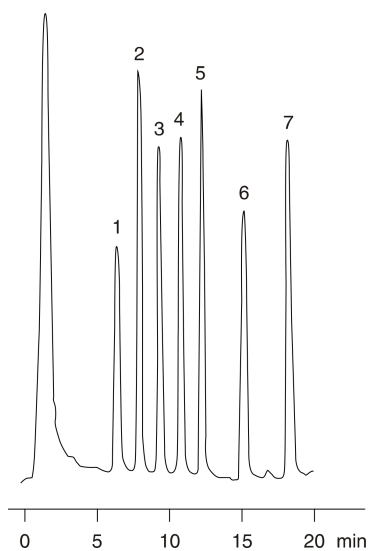

### ANALYTES:

- 1 Ethylamine
- 2 Propylamine
- 3 Butylamine
- 4 Pentylamine
- 5 Hexylamine
- 6 Octylamine
- 7 Decylamine

